

PRIDE[®]
MOBILITY
LIVE YOUR BEST[®]

POWER CHAIR COLLECTION





3
Jazzy® 613 EVO Series

4
Jazzy Air® 2

7
Jazzy® Passport

8-9
Go Chair® / Jazzy® 614 HD

10-11
Jazzy® Elite 14 / Jazzy® Elite HD

12-13
Accessories / Technology

14-15
Specifications

Up To
21
Miles Per Charge*

22"
Turning Radius

Up To
5.1mph
Maximum Speed


300lbs.
Weight Capacity



* 21 mile range for EVO 613 Li; 17 mile range for EVO 613

Jazzy
EVO 613 SERIES

The Jazzy® EVO 613 series provides the freedom to navigate small spaces and tight corners indoors with a narrow base width of only 22". The Jazzy® EVO 613 Li uses a lithium-ion battery that provides longer life and more convenient charging.

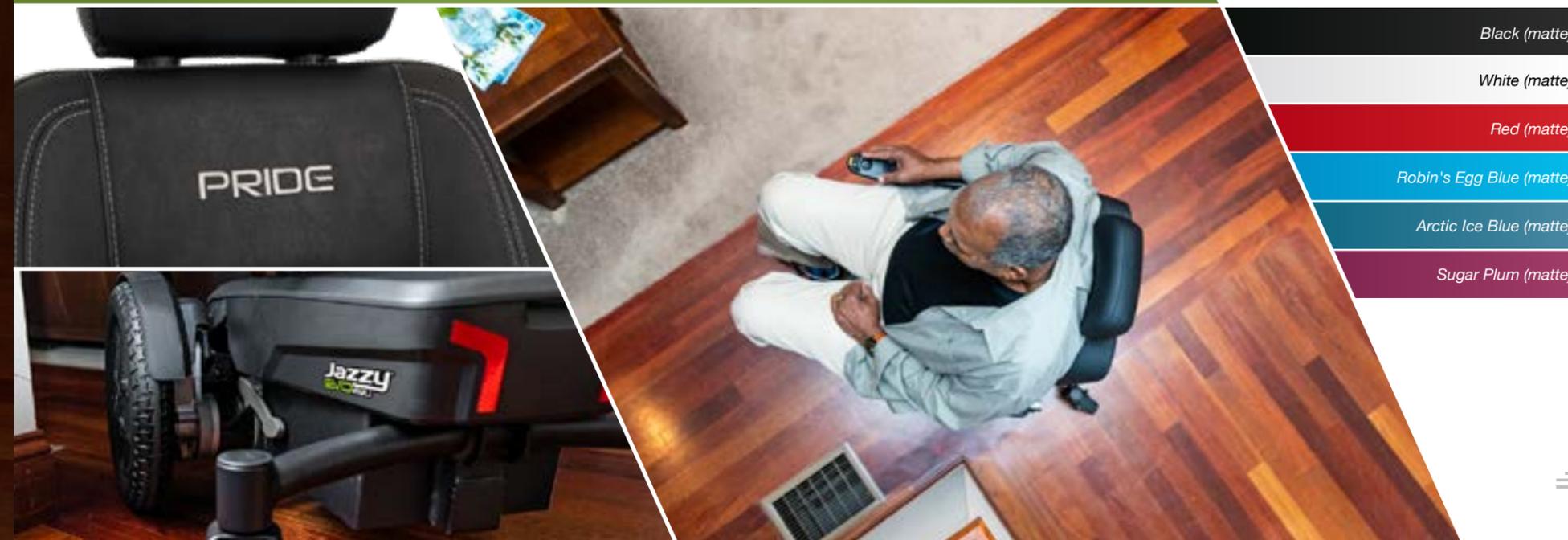

Black (matte)

White (matte)

Red (matte)

Robin's Egg Blue (matte)

Arctic Ice Blue (matte)

Sugar Plum (matte)

Jazzy air₂

Press, elevate and go!
The Jazzy Air[®] 2 elevates 12 inches at the touch of a button. Feel empowered with this sleek technology that allows you to be face to face with your loved ones. Drive at 4 mph while elevated. It's the future of mobility, and the future is yours.



Up To
16.6
Miles Per Charge

20.75"
Turning Radius

Up To
4mph
Maximum Speed
Elevated and Non-Elevated

300lbs.
Weight Capacity

Silver

Pearl White

Ruby Red

Onyx Black

Pink Topaz

Sapphire Blue

Citrine Yellow

Garnet Red
(matte)

Tanzanite
(matte)

Black Pearl
(matte)





Up To
8.6
Miles Per Charge at 220 lbs.

24.25"
Turning Radius

Up To
3.6mph
Maximum Speed

250lbs.
Weight Capacity



Jazzy[®] PASSPORT

Your destination awaits! Whether you are an avid jet setter to remote areas, or a train commuter to the big city, the Jazzy[®] Passport folds and unfolds in just a few steps. Compact and lightweight for hassle-free travel, enjoy abundant storage and comfortable seating.



Black

GOCHAIR®

Compact. Bold. Endless possibilities. The Go Chair® is #1 for portability for easy disassembly in your fast-paced world. Enjoy under seat storage and complete comfort. Choose from six gorgeous colors and let your personality shine through.



Up To
13.2
Miles Per Charge at 200 lbs.

25.5"
Turning Radius

Up To
3.7mph
Maximum Speed

300lbs.
Weight Capacity

Up To
14
Miles Per Charge

20"
Turning Radius

Up To
4.25mph
Maximum Speed

450lbs.
Weight Capacity



Jazzy® 614 HD

Power and precision combine into one exceptional power chair! The Jazzy® 614 HD features outstanding suspension and a tight turning radius. Large drive wheels and six-inch casters tackle uneven terrain, no matter where your day takes you.

Ruby Red

Pearl White

Amber Orange

Citrine Yellow

Rose Quartz

Sapphire Blue

Candy Apple Red

Viper Blue

Jazzy® ELITE 14

An elite design delivers elite performance. The Jazzy® Elite 14 features front-wheel drive for the tightest turns. Large drive wheels provide superior climbing power over life's obstacles. A high-back seat ensures absolute comfort all day.



Up To
18
Miles Per Charge

24"
Turning Radius

Up To
4mph

300lbs.
Weight Capacity

Up To
14.4
Miles Per Charge

24"
Turning Radius

Up To
4mph
Maximum Speed

450lbs.
Weight Capacity



Jazzy® ELITE HD


Experience the power and versatility of the Jazzy® Elite HD. Large front wheels absorb and climb, easily tackling life's everyday obstacles. An intuitive hand control ensures effortless operation and simple maneuvering.

Jazzy Red

Jazzy Blue

Jazzy Red

Jazzy Blue

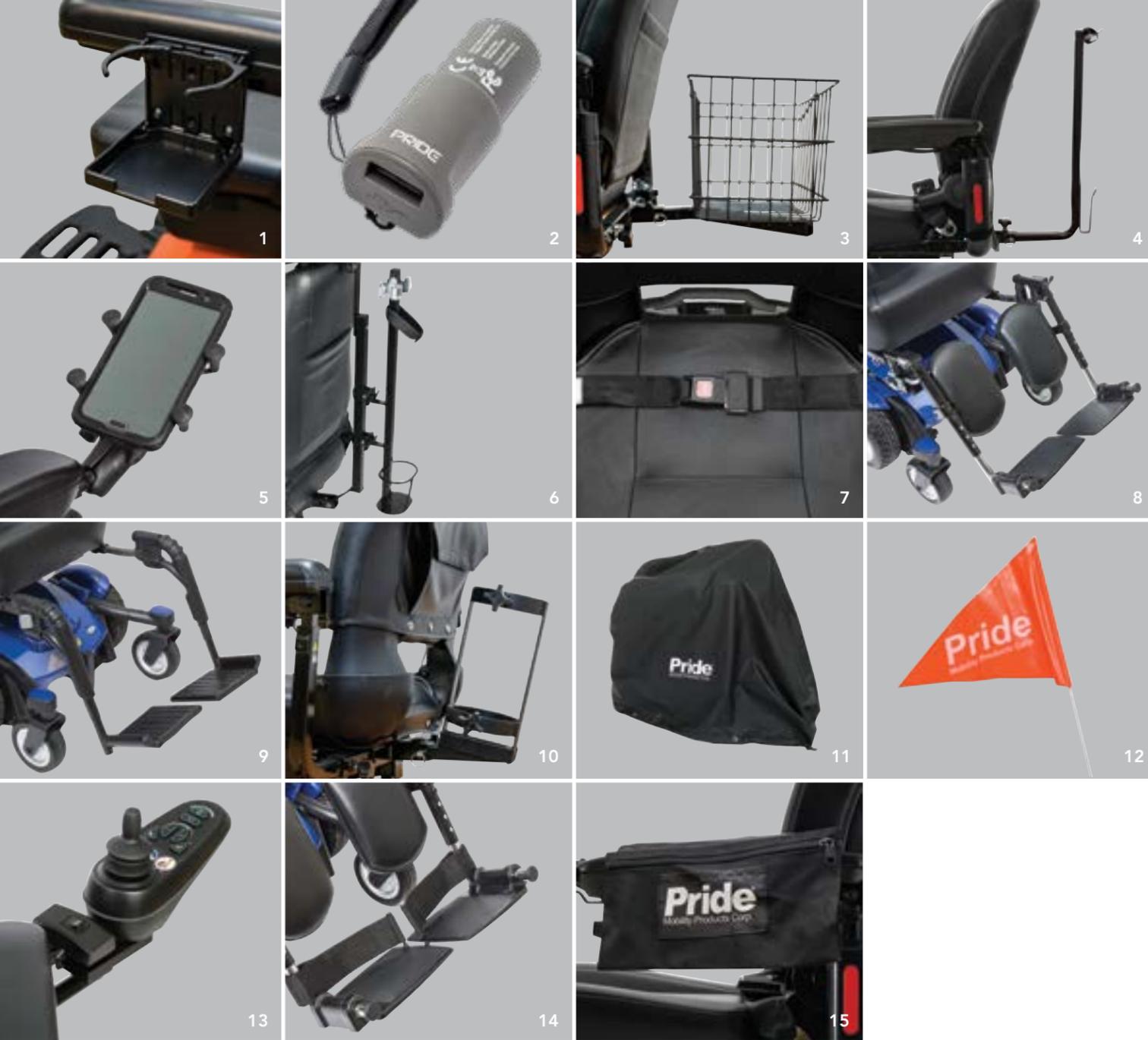

ACCESSORIES

We take innovation to the next level with a wide variety of technologies and state-of-the-art development.

Mid-Wheel 6® Technology and Active-Trac® ATX Suspension allow the active front casters to lift when climbing a threshold or obstacle, while also dropping down for stability when descending.

In-line motors deliver enhanced torque, battery range and reliability for the most dependable performance.

Power elevated seating on the Jazzy Air® 2 features Extra Stability Technology®, mechanically stabilizing the power base, for ultimate security when driving at 4 mph while elevated.

INNOVATIVE TECHNOLOGY



1. Cup Holder | 2. XLR USB Charger | 3. Rear Basket | 4. Walker Holder | 5. Cell Phone Holder | 6. Cane/Crutch Holder
 7. Lap Belt | 8. Elevating Legrests | 9. Swing-Away Legrests | 10. Oxygen Tank Holder | 11. Weather Cover | 12. Safety Flag
 13. Swing-Away Joystick | 14. Heel Loops | 15. Saddle Bag

MODEL NAME	JAZZY® EVO 613 Li	JAZZY® EVO 613	JAZZY AIR® 2	JAZZY® PASSPORT
WEIGHT CAPACITY	300 lbs.	300 lbs.	300 lbs.	250 lbs.
TURNING RADIUS ³	22"	22"	20.75"	24.25"
OVERALL LENGTH /WIDTH ³	35" / 22" (without foot rigging)	35" / 22" (without foot rigging)	41.93" / 23.19"	40" / 23.5"
MAXIMUM SPEED ¹	Up to 5.1 mph	Up to 5.1 mph	Up to 4 mph (elevated and non-elevated)	Up to 3.6 mph
GROUND CLEARANCE ³	3"	3"	1"	3"
FRONT WHEELS	5" solid	5" solid	6" solid	8" solid casters
DRIVE WHEELS	13" solid	13" solid	10" solid	11.75" Perforated
REAR WHEELS	5" solid	5" solid	6" solid	N/A
DRIVETRAIN	Two-motor, in-line, Mid-Wheel 6®	Two-motor, in-line, Mid-Wheel 6®	Two-motor, Mid-Wheel 6®	Two-motor, rear-wheel drive
BRAKING SYSTEM	Regenerative and electro-mechanical	Regenerative and electro-mechanical	Regenerative and electro-mechanical	Regenerative and electro-mechanical
SUSPENSION	Active-Trac® (ATX)	Active-Trac® (ATX)	Active-Trac® (ATX) with Mid-Wheel 6® Technology	Limited
HIGH-BACK SEAT (MAX. SIZE)	20"W x 20"D	20"W x 20"D	20"W x 22"D	N/A
HIGH-BACK SEAT-TO-FLOOR RANGE	High-back: 19.75"-21.25"	High-back: 19.75"-21.25"	21.25-33"	N/A
SYNERGY® SEATING	N/A	20"W x 20"D	N/A	N/A
SYNERGY® SEAT SEAT-TO-FLOOR RANGE	N/A	Synergy : 17.5"-18.5"	N/A	N/A
POWER ELEVATING SEAT ⁶	N/A	N/A	Standard	N/A
STANDARD ELECTRONICS	40A, PG GC3	40A, PG GC3	40A, PG VR2	50A, PG VSI
BATTERY SIZE ⁴	(1) lithium-ion phosphate 25.6V, 32.4Ah	(2) U1 batteries	(2) 12 volt, U-1 (std.), 40 AH (optional)	(1) 25.2V 18 AH lithium ion, (1) 25.2V 10.4 AH lithium ion travel battery
STANDARD BATTERY CHARGER	Off-board, 3.5 amp	Off-board, 3.5 amp	Off-board, 5 amp	Off-board, 2 amp
RANGE PER CHARGE ^{1,2}	Up to 21 miles at 200 lbs. Up to 17.6 miles at 300 lbs.	Up to 17 miles at 200 lbs. Up to 14.4 miles at 300 lbs.	Up to 16.6 miles at 220 lbs. Up to 15 miles at 300 lbs.	18 Ah: Up to 8.6 miles at 220 lbs. 18 Ah: Up to 7.9 miles at 250 lbs. 10.4 Ah: Up to 5 miles at 220 lbs. 10.4 Ah: Up to 4.6 miles at 250 lbs.
BATTERY WEIGHT ⁴	18.4 lbs. total	22 lbs. each	24.5 lbs. each (U-1) 32.5 lbs. each (40 AH)	7.6 lbs., 5.2 lbs. (travel battery)
WEIGHT OF HEAVIEST PIECE ³	94 lbs.	94 lbs.	202 lbs. (with seat)	56.6 lbs. (unit doesn't disassemble) ⁷
STANDARD SEAT WEIGHT	49.8 lbs.	49.8 lbs.	41 lbs.	N/A
STANDARD SEAT SIZE	16" x 16-18", 18" x 18-20", 20" x 18-20"	16" x 16-18", 18" x 18-20", 20" x 18-20"	16" x 16-18", 18" x 18-20", 20" x 18-20"	17.5" x 16.5"
WARRANTY	Frame: 5-year limited Drivetrain: 24-month limited Electronics: 24-month limited Lithium battery: 3-year limited	Frame: 5-year limited Drivetrain: 24-month limited Electronics: 24-month limited SLA battery: 13-month limited	Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months	Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months

GO CHAIR®	JAZZY® 614 HD	JAZZY® ELITE 14	JAZZY® ELITE HD
300 lbs.	450 lbs.	300 lbs.	450 lbs.
25.5"	20"	24"	24"
26" / 22.5" (without foot rigging)	35.5" / 25.25"	35.5" / 24.25"	35.5" / 24.25"
Up to 3.7 mph	Up to 4.25 mph	Up to 4 mph	Up to 4 mph
1.75"	1.625"	3"	3"
6" solid casters	6" solid casters	3" solid anti-tips	3" solid anti-tips
9" solid	14" solid	14" solid	14" solid
3" solid anti-tips	6" solid casters	6" solid casters	6" solid casters
Two-motor, rear-wheel drive	Two-motor, in-line, Mid-Wheel 6®	Two-motor, in-line, front-wheel drive	Two-motor, in-line, front-wheel drive
Regenerative and electro-mechanical	Regenerative and electro-mechanical	Regenerative and electro-mechanical	Regenerative and electro-mechanical
Limited	Active-Trac® (ATX)	Limited	Limited
N/A	24" W x 24" D	22"W x 22"D	24" W x 24" D
N/A	20.5-23.5"	20-22"	20-22"
N/A	24" W x 24" D	20"W x 20"D	24"W x 24"D
N/A	17.75-19.75"	15.75-17.75"	15.75-17.75"
N/A	N/A	N/A	N/A
40A, PG GC3	40A, Dynamic LiNX	40A, Dynamic LiNX	40A, Dynamic LiNX
(2) 18 AH	(2) NF-22	(2) NF-22	(2) NF-22
Off-board, 2 amp	Off-board, 5 amp	Off-board, 5 amp	Off-board, 5 amp
Up to 13.2 miles at 200 lbs. Up to 8.7 miles at 300 lbs.	Up to 14 miles	Up to 18 miles	Up to 14.4 miles
31 lb. battery box	37 lbs. each	37 lbs. each	37 lbs. each
36 lbs.	122 lbs.	99 lbs.	99 lbs.
30 lbs. (compact foldable)	52 lbs. (high-back)	52 lbs. (high-back)	52 lbs. (high-back)
18" x 17"	20" x 18-20", 22" x 20-22"	18" x 18-20", 20" x 18-20", 22" x 20-22"	20" x 18-20", 22" x 20-22"
Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months	Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months	Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months	Frame: 5-year limited Drive motors: 24 months Electronics: 24 months Batteries: 13 months

SPECIFICATIONS

PRIDE[®]

MOBILITY
LIVE YOUR BEST[®]

Beamsville, ON | 888-570-1113 | pridemobility.ca |    



Note: Transport of the Jazzy[®] Passport item on commercial aircraft is subject to federal regulations and airline policy. Please consult your airline prior to your flight to ensure your Jazzy[®] Passport and its lithium-ion battery are permitted.

- 1 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition, and tire condition. This specification can be subject to a variance of (+ or -) 10%.
- 2 Tested in accordance with ANSI/RESNA, MC Vol2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- 3 Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or -) 3%.
- 4 Battery weight may vary based on battery manufacturer.
- 5 AGM or gel cell type required.
- 6 Maximum weight capacity is reduced to 400 lbs. (181 kg) on 18" deep seat configurations.
- 7 Seat weight included. Unit does not disassemble.
- 8 40A, PG 6-key VR2 required with power elevating seat option.
- 9 Not all accessories are available for all products. Please speak with your Provider for details and availability.
- 10 High-back only (STF height: 22.5-26.5"); Range per charge: Up to 13 miles; Base weight: 114.5 lbs.

NOTE: All specifications subject to change without notice.

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Dimensions listed are for power chair only. Overall measurements will vary based on seating and accessory selections.