

Model	Maxima
Part Number	3W: SC901 4W: SC941
Weight Capacity	500 lbs.
Maximum Speed¹	Up to 5.8 mph
Ground Clearance²	3W: 2.5" at motor; 4W: 2.25" at motor
Turning Radius²	3W: 49"; 4W: 53"
Overall Length²	47.5"
Overall Width²	27" (not including armrest)
Seat-to-Ground Height²	(3W/4W): 22.1" – 24.1" Power Elevating Seat: 22.32" – 28.57"
Seat-to-Deck Height²	(3W/4W): 17.13" – 19.13" Power Elevating Seat: 17.36" – 23.61"
Front Tires²	3W: 10" x 3.6" solid 4W: 10" x 3.6" solid
Rear Tires²	10.75" x 3.6" solid
Suspension	N/A
Range Per Charge^{1,4}	3W: Up to 19.8 miles 4W: Up to 19.5 miles
Total Weight With and Without Batteries⁶	3W (w/out batteries): 183.5 lbs. 3W (w/batteries): 256.5 lbs. 4W (w/out batteries): 196.8 lbs. 4W (w/batteries): 269.8 lbs.
Heaviest Piece When Disassembled²	Unit (w/out seat and batteries) 3W: 130.3 lbs. 4W: 143.6 lbs.
Standard Seating²	Type: Depth-adjustable, high-back Weight: 53.2 lbs. Material: Black vinyl Dimensions: 22" x 18"-20" (depth adjustable) Back height: 19.75"
Drive System	Sealed transaxle, 24-volt DC motor
Dual Braking System	Regenerative and electromechanical
Battery Requirements^{3,5}	Size: (2) Gel, 12V, 58Ah, NF-22 Weight: 36.5 lbs. each
Battery Charger	Off-board, 5-amp
Warranty	Frame: lifetime limited Drivetrain: 2-year limited Electronics: 2-year limited Batteries: 13-month limited
Color Option	Garnet Red (matte) Ocean Blue (matte)

¹ Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).

² Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³ AGM or gel cell type required.

⁴ Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

⁵ Battery weight may vary based on manufacturer.

⁶ Includes standard seat weight.

PRIDE[®]
MOBILITY
LIVE YOUR BEST[®]

Beamsville, Ontario • 888-570-1113

pridemobility.ca |    

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

0002977-05-02-22

MAXIMA

ULTRA HEAVY-DUTY FOR MAXIMUM PERFORMANCE



LIVE YOUR BEST[®]
pridemobility.ca

HEAVY-DUTY PERFORMANCE



Full LED light package

Infinite position tiller lever (mech lock)



ALL NEW REDESIGN!

NEW! Tiller adjustment

- Infinite position tiller lever (mech lock)

NEW! Seat

- 22" x 18"-20" High back depth adjustable
- 106° limited recline
- 4-point swivel (no sliders)
- 7" headrest with 3" of adjustment

NEW! USB charger built into the tiller for convenient charging of smart phones and portable electronic devices

NEW! Full LED light package includes front headlights, turn signals and rear taillights

NEW! Ambient voltmeter screen brightness automatically adjusts based on light intensity


NEW! Delta tiller with ergonomic wraparound handles lets you operate the scooter with one hand

- Longer per-charge range
- Standard with front and side bumper (3W) and side bumper (4W)
- XLR battery charger on tiller

Console

- Ambient voltmeter
- On/off key switch
- Horn
- Speed control
- Light switch
- Turn signal button

Colors

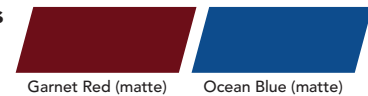

Options

- Power elevating seat available on 3 & 4 wheel
- 6.25" of additional height to access cabinets, countertops and other out-of-reach items

Accessories

- Additional rearview mirror
- Rear bumper (3w, 4w)
- Cup holder
- Saddle bag
- Forearm crutch holder
- Safety flag
- Front bumper (4w)
- Single/double cane/crutch holder
- Oxygen holder
- Weather covers
- Quad cane holder
- Wishbone crutch holder
- Rear basket