

Jazzy® Carbon¹

| | |
|----------------------------------|---|
| HCPCS Code | K0899 Power Mobility Device, not coded by DME PDAC or does not meet criteria |
| Weight Capacity | 300 lbs. |
| Turning Radius ³ | 24" |
| Overall Width ³ | 24" |
| Overall Length ³ | 37" |
| Maximum Speed ² | Up to 3.7 mph |
| Ground Clearance ³ | 2" |
| Front Casters ³ | 7" |
| Drive Wheels ³ | 8.5" solid |
| Rear Casters ³ | N/A |
| Drivetrain | 24v/250W brushless motors |
| Braking System | Regenerative and electromechanical |
| Suspension Type | Front |
| Standard Seat ³ | 17" x 17" black foam |
| Seat-to-Floor Range ³ | 20.25" with standard cushion 18.75" without cushion |
| Per Charge Range ^{3,5} | Up to 9.3 miles |
| Standard Electronics | 35A Brushless Micon |
| Component Weights ⁶ | Base: 39.6 lbs. Battery: 4 lbs. |
| Total Weight ⁴ | 43.6 lbs. with batteries |
| Battery Requirements | 24V 12AH (288WH) lithium-ion |
| Battery Charger | 2 amp, lithium-ion off-board |
| Warranty | Lifetime limited on frame; 1-year limited on electrical control/lever system, motor/drive system, bearing sleeve and battery |



FDA Class II Medical Device¹

PRIDE[®]
MOBILITY

800-800-8586 | pridemobility.com



HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

102124

- Jazzy® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.
- Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition. These specifications can be subject to a variance of (+/- 10%).
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- Total weight includes base, seat and batteries.
- Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- Battery weight may vary based on manufacturer.

Jazzy[®] carbon



**FOLDS FLAT
FOR TRANSPORT**

**OPENS
WITH
ONE-HAND
SNAP
ACTION!**



**LIGHT
ENOUGH TO
PICK UP!**

AMERICA'S #1 TRAVEL POWER CHAIR

pridemobility.com |    

FDA Class II Medical Device¹

Jazzy[®] carbon

39.6
LBS
WEIGHT WITHOUT BATTERY

Padded seat back with storage pocket

Quick disconnect joystick with integrated USB charger

Mesh Cup holder


Ultra-lightweight, carbon fiber frame

Compact base easily folds and unfolds for easy transport in a few easy steps

Airline compliant lithium battery*

Front suspension for a smooth ride

Under-seat storage

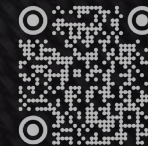

FDA Class II
Medical Device¹

UNFOLD, SNAP, AND GO!

This is our most portable Jazzy[®] yet. In just seconds, the Jazzy Carbon folds flat for easy transport, and unfolds again to get you moving. Ultra-light but exceptionally strong, this power wheelchair is your ultimate travel companion.

Weighing less than 40 lbs. (without battery), the Jazzy Carbon still provides top speeds up to 3.7 mph with a weight capacity of 300 pounds.

FEATURES

- Quick disconnect joystick with built in USB charger
- 24V 12Ah (288WH) airline compliant lithium battery*
- Front suspension for a smooth ride
- Ultra-lightweight, carbon fiber frame specifically designed for traveling
- Compact base easily folds and unfolds for transport in a few simple steps
- Padded seat back with storage pocket and soft, removable seat cushion for all-day comfort
- Free-wheel levers easily disengage the brake, allowing this lightweight power chair to be pushed like a manual chair
- Under-seat storage and adjustable lap belt
- Standard mesh cup holder

AVAILABLE COLORS



Black

White

ACCESSORIES

- Additional battery
- Attendant control bracket
- **NEW!** Battery bag
- Headrest
- **NEW!** Jazzy[®] Travel Kit (includes a bag to house the folded unit and a separate battery bag)
- Military patches
- Stealth[®] Simplicity cushion

NEW!



Battery Bag



Height: 28.3" | Width: 24" | Depth: 12.2"



The Jazzy[®] Travel Kit includes a bag to protect your folded Jazzy[®] Carbon and a separate bag for the battery.

Folds quickly and easily for transport or storage.

Front suspension provides a smooth ride.

* Transport of this item on commercial aircraft is subject to federal regulations and airline policy. Please consult your airline prior to your flight to ensure your Jazzy[®] Carbon and its lithium ion battery are permitted.