

# SSAPT B4

Simple Solutions Adapter Interface




# SSAPT B4

## Simple Solutions Adapter Interface

### OVERVIEW

The Simple Solutions Adapter Interface (**SSAPT B4**) is designed to be used on Quantum's TB-series of chairs.

Stealth's goal was to design a mounting system that would give the ability to mount a wide variety of back options with one key feature in mind, simple.

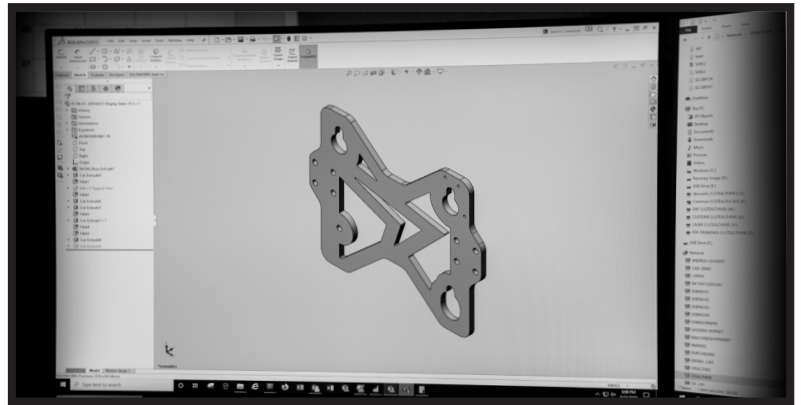

#### **CENTRAL MOUNTING** 1

The center Plate mounts directly to the TB-series center spine.

#### **WING WIDTH & ADJUSTABILITY** 2

Available 3 7/8" in width adjustment per wing, allowing for a wider range of backrest widths and lateral offset.

#### **MOUNTING POINT ADJUSTABILITY** 3

360 degrees of rotation is achievable by using increments of 90 degrees.

#### **QUICK RELEASE LEVER** 4

Locks the SSAPT B4 in place & allows easy, quick removal of the back when the thumb screw is not being utilized.


#### **STEALTH UniLink™ INTEGRATION** 5

The SSAPT B4 comes with 4 UniLink™ sections.


#### **UniLink™ ADJUSTABILITY** 6

Available 3" of height adjustment per UniLink™ section.

- Adjustment of the links can effect:
- Depth
  - Back Rotation
  - Back Angle
  - Back Height

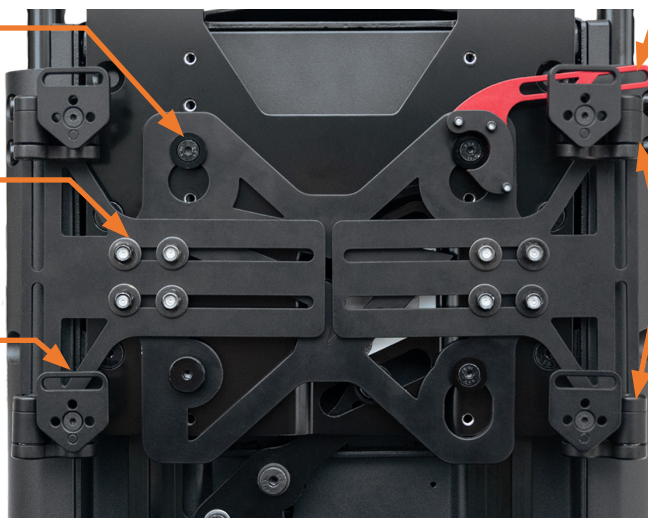

### FEATURES

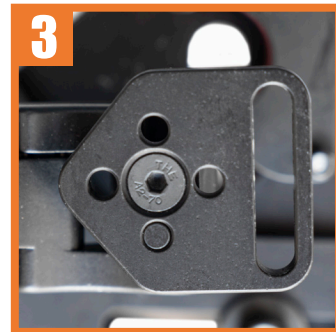

The SSAPT<sub>B4</sub> was designed to work with the following brands:

Stealth Products®

Invacare®

Freedom Designs

Symmetric Designs

Comfort Company

Matrx® Backs

Jay

Roho®

Ride®

Axiom Custom Seating

AES® Custom Seating

**COMPATIBILITY**



SSAPT<sub>B4</sub> shown with Stealth Tarta Original Custom Back

## IMPORTANT TO NOTE

Contact a Stealth CSR if a specific brand is not listed. An unlisted back may still be compatible.

In standard configuration, the SSAPT<sub>B4</sub> is compatible with backs that have mounting points between 9¾ and 20⅝ inches.

## FOR MORE INFORMATION

### SCAN THE CODE



[WWW.STEALTHPRODUCTS.COM](http://WWW.STEALTHPRODUCTS.COM)

**OR**



[Productorders@pridemobility.com](mailto:Productorders@pridemobility.com)

+1 (888) 570-1113

Distributed by  
Pride® Mobility Products Company  
5096 South Service Road  
Beamsville, ON L0R 1B3



**Stealth Products**  
*Solutions For Life*®

CASSAPT4SS0820