



### TRU-Comfort 2 Back and Seat Cushion

- Adjustable positioning and pressure management
- Compatible with both TRU-Balance® 3 Power Positioning Systems and the Static Seat Frame System
- Available with four-way stretch, recovery, fluid-resistant, breathable cover
- Additional 1" of viscoelastic foam for maximum comfort and pressure management



### TRU-Balance® 3 Static Seating (for Edge 2.0 only)

- Easily adjustable and versatile design that can grow with the client
- Easy customization
- 300 lb. weight capacity
- Accepts a wide variety of after-market seats and backs
- 20° of limited recline adjustment (-5 / +15) via gas strut within back cover
- Available in 12 back shroud colors



### TRU-Balance® 3 Manual Recline (for Edge 2.0 only)

- Solid curved back pan and 2.5" raised pivot to reduce shear forces
- Attendant activated recline seat
- Width: 12" - 24"; depth: 12" - 24"
- 300 lb. weight capacity
- 175° of recline
- Can be ordered with TRU-Balance® 3 Tilt and TRU-Balance® 3 Static Seating
- Available in 12 back shroud colors

### Complimentary Accessory



Cup Holder  
(front mounted)

### Premium Accessories (Ask a Quantum Rehab® Specialist for a complete list.)



USB Port

Glove Box

Clothing Guard

Push Bar

Hydration System

Transfer Bar

# TRU BALANCE<sup>3</sup>

## POWER POSITIONING SYSTEMS

Tailor-Made with  
Conforming Versatility



# TRU BALANCE<sup>3</sup>

## POWER POSITIONING SYSTEMS

- Make fine adjustments while client is in the chair
- Order your seat width independent of your back width at no additional charge
- Internally or externally rotate cantilever armrests for optimal positioning
- Easily retrofit power tilt, recline, and scissor lift
- Easy back height and depth adjustment
- Easy center of gravity adjustment
- Available with a complimentary cup holder
- Premium accessories available
- Tool kit included for quick adjustments



### TRU-Balance<sup>3</sup> Power Tilt

- 50° of tilt
- Lower standard seat-to-floor height
- Easy, unmatched adjustment and durability
- Maintains center of gravity shift
- Available in 12 back shroud colors



### TRU-Balance<sup>3</sup> Power Adjustable Seat Lift with iLevel<sup>®</sup>

- Offers 10" of lift in 16 seconds at up to 3.5 mph with iLevel<sup>®</sup>
- Allows easy access to batteries
- 300 lbs. weight capacity
- Available in 12 back shroud colors



### TRU-Balance<sup>3</sup> Power Recline

- 85° - 175° of recline
- 300 lbs. weight capacity
- Does not increase seat-to-floor height
- Standard shear reduction and cantilever armrests
- Available with TRU-Balance<sup>3</sup> Power Tilt and TRU-Balance 3 Lift
- Available in 12 back shroud colors



### TRU-Balance<sup>3</sup> 10" HD Power Adjustable Seat Lift

- Offers 10" of lift
- 400 lbs. weight capacity
- Seating sizes up to 24" x 25"
- Adjustable seat-to-floor height aids in daily transfers
- Allows easy access to batteries



### TRU-Balance<sup>3</sup> HD Power Tilt

- 450 lbs. weight capacity
- 50° of tilt
- Easy, unmatched adjustment and durability
- Maintains center of gravity shift
- Available in 12 back shroud colors



### TRU-Balance<sup>3</sup> HD Power Recline

- 450 lbs. weight capacity
- 85° - 175° of recline
- Does not increase seat-to-floor height
- Standard shear reduction and cantilever armrests
- Available in 10 back shroud colors
- Available with TRU-Balance<sup>3</sup> HD Power Tilt




### Q-Logic EX Attendant Control

- Full proportional joystick
- Standard on power tilt and power recline when ordered with an Enhanced Display or SCIM Module
- Allows caregivers the ability to provide patient assistance



### Recline Back Interface for After-Market Backs

- Optional integration system allows for an unprecedented level of compatibility between TRU-Balance<sup>3</sup> Power Recline and most after-market back rests
- Retains full functionality, including shear reduction

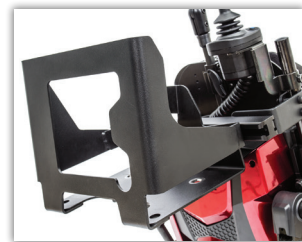

### Vent Tray Solutions

- Numerous vent tray and oxygen holder configurations are available like our Trilogy Vent Tray (for Trilogy 100 or 200 ventilators) and our compact vent tray (for LTV style ventilators)
- The Quantum Rehab<sup>®</sup> Applications Team is available for many specialized vent tray modifications. Contact the Rehab Applications Team for more information.



### TRU-Balance<sup>3</sup> Power Articulating Foot Platform

- Broad articulation range for enhanced adjustment
- 450 lbs. weight capacity
- Unmatched adjustment, durability and strength
- Multiple calf pads, calf panel and footplate options



### TRU-Balance<sup>3</sup> Power Articulating Leg Rests

- 70° angle
- Actuator is mounted to the outside of leg rest without increasing seat-to-floor height
- Connected with a ¼" headphone-style jack connector (pull out to release)
- Depth and height-adjustable calf pads
- 4" or 6" articulation
- Swing-away